


BURETTES - "schilling" - "class AS"

- Described as automatic burettes with fast titration and practical usage.
- Especially suitable for field analysis where there is a risk of breakage.
- Available 2 types:
 - clear burette with blue scale
 - amber burette with white scale
- Calibrated to deliver (Ex) and the tolerances are according to ISO 385 Class "AS" standards.
- The entire system consists of a burette tubing with automatic zeroing system together with a valve stopcock, a reservoir bottle and a bottle stand.
- Rapid filling by squeezing the bottle and fast titration through valve stopcock.
- Inscriptions and graduations are printed either with blue enamel or with white enamel which is specially designed for glass.
- Clear schilling burettes are manufactured with schellbach stripe. Amber schilling burettes are manufactured without schellbach stripe.
- 10 ml burette has a reservoir bottle of 500 ml where as 25 and 50 ml burettes are used with bottles of 1000 ml.
- Other technical specifications are exactly the same as other types of burettes.

clear with schellbach	amber w/o schellbach	volume	subdivision	pack quantity
020.03.010	020.05.010	10 ml	0.05 ml	1 piece
020.03.025	020.05.025	25 ml	0.10 ml	1 piece
020.03.050	020.05.050	50 ml	0.10 ml	1 piece

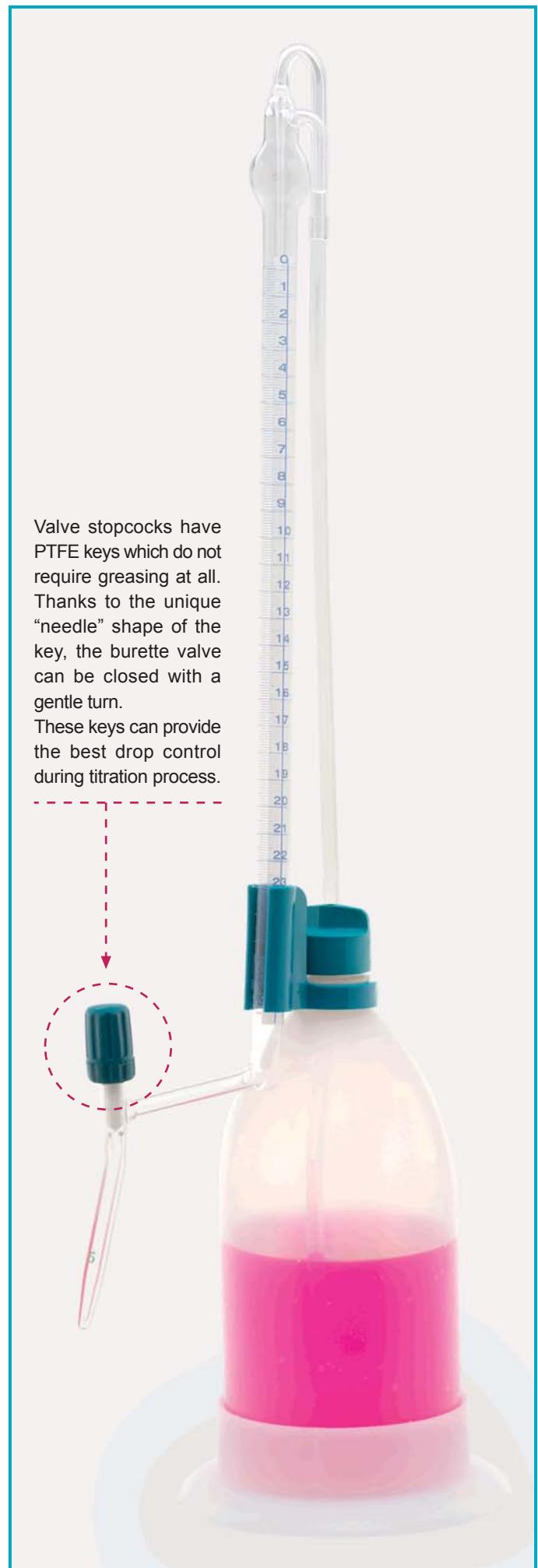

Valve stopcocks have PTFE keys which do not require greasing at all. Thanks to the unique "needle" shape of the key, the burette valve can be closed with a gentle turn. These keys can provide the best drop control during titration process.