

BLOOD LANCETS - "manual"

- Manufactured from medical grade stainless steel.
- The gamma-sterilized lancets are individually packed.

Lancets are very incisive and have very sharp edges which penetrates into the skin very easily.



catalogue number	quantity in pack	quantity in box
088.01.001	200 pieces	50 pack

BLOOD LANCETS - "automatic"

- These lancets are user safe since the needle is fully shielded before and after use.
- They have a pre-loaded, button operated design.
- These lancets are for single-use only, sterilized by gamma irradiation.
- Patient comfort due to high quality, ultra sharp tri-bevel needle, high speed penetration of 1.8 mm penetration depth.




catalogue number	quantity in pack	quantity in box
088.02.003	100 pieces	50 pack



Twist off the protective cap/safety cover



Hold the lancet firmly against selected finger site.

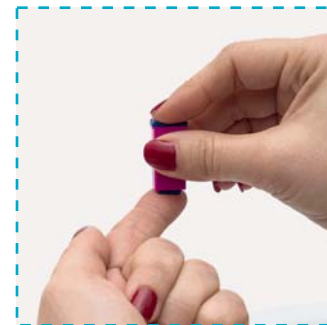

Press the release button to activate the device.



Gently massage the finger to obtain the required blood.